



BIOTECH AB

WWW.BIOTECH-SWEDEN.COM

Application Data

No. 1051A

ChromaNik
ChromaNik Technologies Inc.


エリストマイシンエステルとの分離

Sunniest C18, 5 μm 250 x 4.6 mm i.d.

Erythromycin Estolate

Chromatographic Purity Test

Source Reference USP 30

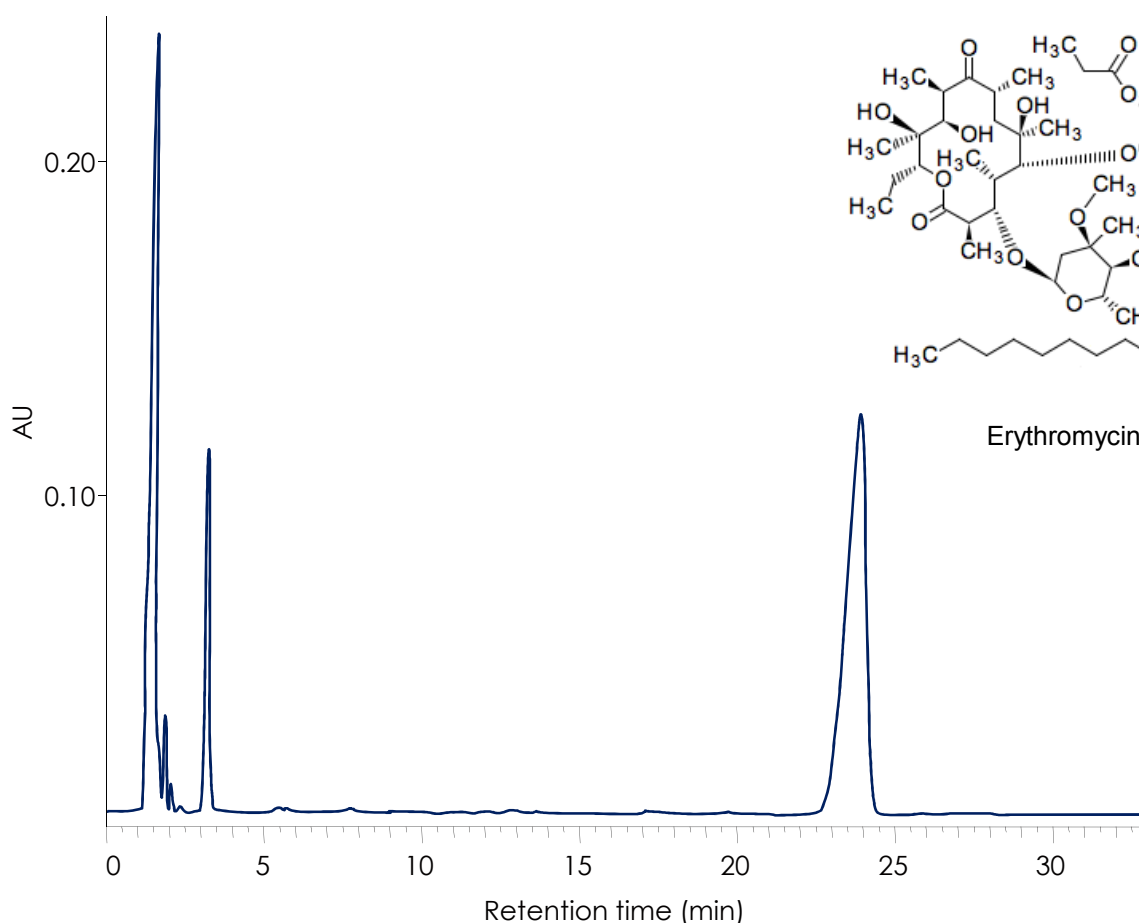

Column: Sunniest C18, 5 μm 250 x 4.6 mm i.d.

Mobile phase: 50ml of 35gm/lit dibasic potassium phosphate, pH 8.0 with dilute phosphoric acid, add 400 ml water add 165 ml 2methyl 2propanol add 30ml ACN, adju volume to 1lit with water

Flow rate : 2.0. mL/min

Detection : UV@215 nm

Column Temperature : 65 °C

Sample Temperature : Ambient

Injection Volume : 200 μ L

Biotech AB in partnership with ChromaNik Technologies

Biotech AB (Sweden) Email: info@biotech.se

Biotech USA Email: info@biotechusa.us

Biotech India Pvt Ltd Email: tech@biotechltd.in

ChromaNik Technologies Inc. (Japan) Email: info@chromanik.co.jp

No. 1051A