

IC Columns Introduction

Alltech® IC Column Specifications						
Column	Composition	pH Range	Applications	EPA Methods	Optimized for Suppressed or Non-Suppressed Conductivity	Page
Anion Exchange Columns						
Allsep™	polymer-based anion exchanger, 7µm	pH 2–10	inorganic anions, weak and strong acid ions, metal complexes, organic acids	300.0, Part A	both	99
Novosep™ A-2	polymer-based anion exchanger, 5µm	pH 3–12	inorganic anions and oxyhalides; ideal for separation of seven common anions plus three oxyhalide anions in one run	300.0, 300.1, 317.0, 326.0	suppressed	100
AN1™	polymer-based anion exchanger, 9µm	pH 2–13	inorganic anions, weak and strong acid ions, organic acids	—	both	100
Anion/S	silica-based anion exchanger, 10µm	pH 2–5.5	inorganic anions	—	non-suppressed	101
Anion/R	polymer-based anion exchanger, 10µm	pH 2–12	inorganic anions, weak and strong acid ions	—	non-suppressed	101
PRP™-X100	polymer-based anion exchanger, 10µm	pH 1–13	inorganic anions, weak and strong acid ions	—	non-suppressed	103
Cation Exchange Columns						
Universal Cation	silica-based cation exchanger, 7µm	pH 2–7	groups I and II cations, amines, divalent transition metals	300.7	both	102
Universal Cation HR	silica-based cation exchanger, 3µm	pH 2–7	groups I and II cations, amines, divalent transition metals; smaller particle size for improved peak resolution	300.7	both	102
PRP™-X200	polymer-based cation exchanger, 10µm	pH 1–13	groups I and II cations (separate runs), amines, lanthanides	—	non-suppressed	103
Organic Acid Columns						
PRP™-X300	polymer-based sulfonated cation exchanger	pH 1–13	simple aliphatic carboxylic acids and alcohols	—	non-suppressed	103
Anion Exclusion	polymer-based cation exchanger, 10µm	pH 1–13	organic acids, weak acid anions	—	non-suppressed	104
Specialty Columns						
Surfactant C8	silica-based reversed-phase, 5µm	pH 2–7	short chain anionic surfactants	—	non-suppressed	105
Surfactant C18	silica-based reversed-phase, 5µm	pH 2–7	long chain or aromatic anionic surfactants	—	non-suppressed	105

IC Guard Kits Cross Reference

Kit includes holder and three guard cartridges.

Guard Reference Table

Brand	Column Type	Guard Cartridge Type	All-Guard™ Kit Part No.
Alltech	Allsep™, Allsep™ A-2	GA-1	38109
	Novosep™ A-2	GA-1	38109
	Anion/R, Anion/S	GA-1	38109
	Universal Cation	GC-2	27109
	Universal Cation HR	GC-3	23115
	Anion Exclusion	GIE-1	38117
	Transition Metal C18	GRP-2	38123
	Surfactant/R	GRP-1	38122
	Surfactant C8, C18	GRP-2	38123
	IonPac® AS4A-SC, AS7, AS11, AS22, AS24, AS26, AS27	GA-1	38109
Dionex	IonPac® CS3, CS10, CS11	GC-1	38113
	IonPac® CS12, CS14, CS15	GC-2	27109
	IonPac® NS1	GRP-1	38122
	IonPac® AS1, AS5	GIE-1	38117

Guard Reference Table (continued)

Brand	Column Type	Guard Cartridge Type	All-Guard™ Kit Part No.
Waters	IC-Pak™-A	GA-1	38109
	IC-Pak™-M/D	GC-2	27109
	Ion Exclusion	GIE-1	38117
Hamilton	PRP™-X100	GA-1	38109
	PRP™-X200	GC-1	38113
	PRP™-X300	GIE-1	38117
Misc.	Sarasep ANIONSEP AN1™	GA-1	38109
	Bio-Rad®, Aminex®, HPX-87	GIE-1	38117
	Interaction™ ION-120	GA-1	38109
	Interaction™ ION-210	GC-1	38113
	Interaction™ ORH-801	GIE-1	38117
	Metrohm® Metrosep	GA-1	38109

more info

Cartridges and holders can also be purchased separately. See specific phase pages for replacement cartridge ordering information.




Alltech® Anion Columns

Allsep™ Anion IC Columns

- Resolves fluoride away from the water dip
- Meets requirements for U.S. EPA Method 300.0 part A

The Allsep™ Anion Column is compatible with common IC mobile phases: carbonate, bicarbonate, p-hydroxybenzoic acid, phthalic acid, succinic acid, and sodium octane sulfonate.

Use this column for performing EPA Method 300.0 Part A "Determination of Inorganic Anions in Water" per U.S. EPA guidelines.

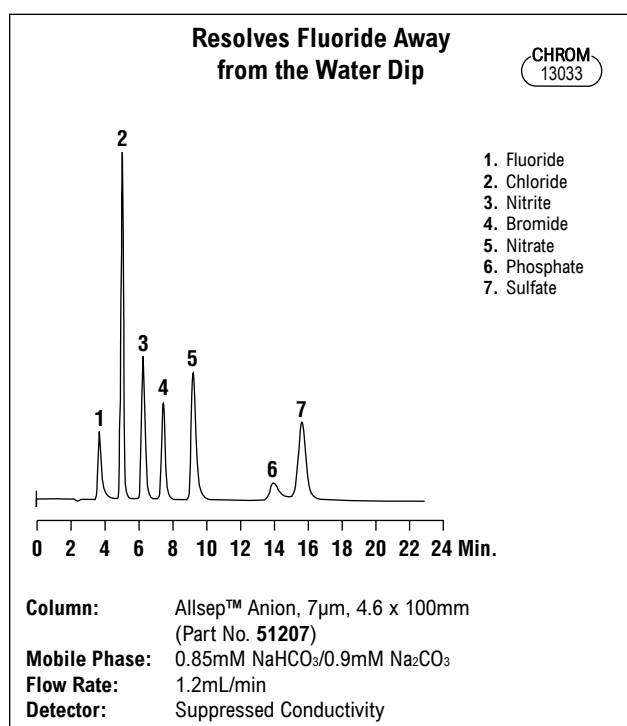

Allsep™ Anion Specifications	
Composition:	Methacrylate based w/quaternary ammonium functional groups
Particle Size:	7µm
Mobile Phase Limits:	pH 2–10, 0–100% organic modifier

Allsep™ Anion Columns			
Packing	i.d. x Length	PEEK Part No.	SS Part No.
Allsep™ Anion	4.6 x 50mm	51213	51214
	4.6 x 100mm	51207	51200
	4.6 x 150mm	51209	51208

related product

Save money with IC All-Guard™ Kits on page 98. Kit includes one holder and three cartridges.

Alltech® and Transgenomic® Anion Columns

Alltech

Novosep™ A-2 Anion IC Columns

- Separates seven inorganic anions and three oxyhalide anions in a single run
- Meets the requirements of U.S. EPA Methods 300.0, 300.1, 317.0, and 326.0



The Novosep™ A-2 Anion Column is ideal for the separation of 10 inorganic anions by suppressor-based ion chromatography using U.S. EPA Method 300.1. The column can also be used to determine inorganic anions in drinking water per U.S. EPA Method 300.0 and trace bromate in drinking water using U.S. EPA Methods 317.0 and 326.0.



Transgenomic® IC Sep AN1™ IC Column

- Separates anions with carbonate/bicarbonate and hydroxide mobile phases

The Transgenomic® IC Sep AN1™ Column is compatible with both suppressor-based and non-suppressed detection. A variety of mobile phases including sodium carbonate/bicarbonate, sodium hydroxide, and borate/gluconate may be used.

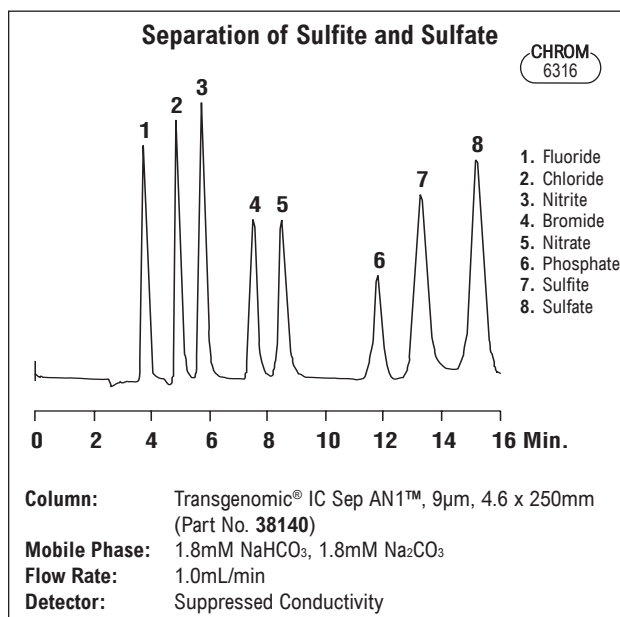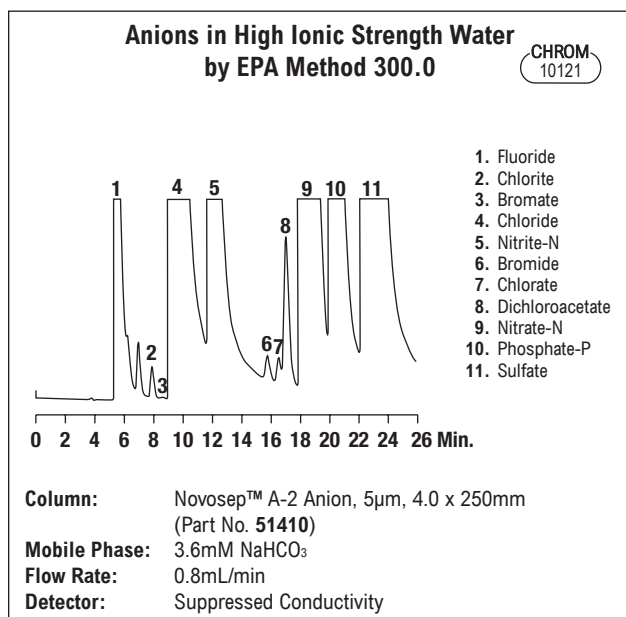

Novosep™ A-2 Anion Specifications

Composition:	Polyvinyl alcohol-based w/quaternary ammonium functional groups
Particle Size:	5µm
Mobile Phase Limits:	pH 3–12, 0–100% organic modifier, 10% modifier is practical upper limit

Transgenomic® IC Sep AN1™ Specifications

Composition:	Poly(styrene-divinylbenzene) alkyl dimethyl ethanol ammonium functional group
Particle Size:	9µm
Exchange Capacity:	0.05meq/g
Mobile Phase Limits:	pH 2–13, No organic solvents

Novosep™ A-2 Anion Column

Packing	i.d. x Length	PEEK Part No.
Novosep™ A-2	4.0 x 250mm	51410

Transgenomic® IC Sep AN1™ Columns

Packing	i.d. x Length	Part No.
AN1™	4.6 x 250mm	38140

Guards for Novosep™ A-2 Anion Columns

Description	Qty.	Part No.
GA-1 Anion Cartridges	3	38108
All-Guard™ Cartridge Holder*	ea	80101

*Direct-Connect™ column coupler included.

Choose a guard column or integral guard discs to protect your IC Sep Column.

Transgenomic® IC Sep Guards

Description	Qty.	Part No.
Guard Column, 50 x 4.6mm	ea	38145
Guard Column Guard-Disc™	5	38146

more applications

To view our complete searchable chromatogram database visit www.discoverysciences.com/chromdb/




Alltech® Anion Columns

Anion/S IC Columns

- Silica-based for symmetrical peak shapes
- Separates inorganic and organic anions

The Anion/S Columns are best suited for routine separations of chloride, bromide, nitrate, and sulfate. Not suitable for fluoride analysis.

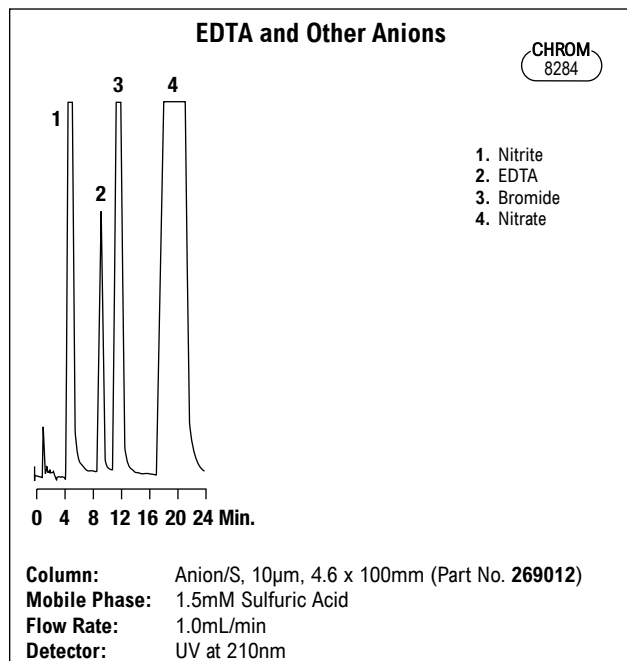

Anion/S Specifications

Composition:	Silica w/quaternary ammonium ion exchanger
Particle Size:	10µm
Exchange Capacity:	0.25meq/g
Mobile Phase Limits:	pH 2–5.5, 0–100% organic modifier

Anion/S Columns

Packing	i.d. x Length	PEEK Part No.	SS Part No.
Anion/S	4.6 x 100mm	269012	269013
	4.6 x 250mm	—	269001


Alltech



Anion/R IC Columns

- Polymer-based for broad pH stability
- Separates seven common inorganic anions

Because of their high pH stability, these columns are ideal for the separation of anions using high pH mobile phases such as p-hydroxybenzoate, sodium hydroxide, and sodium carbonate/bicarbonate. The larger dimension columns (150mm and 250mm) are suitable for a high-resolution separation of the seven common inorganic anions.

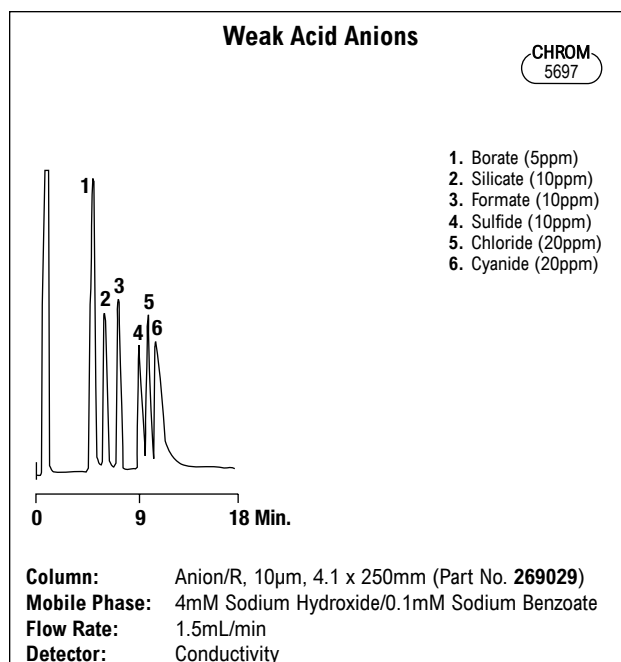

Anion/R Specifications

Composition:	Poly(styrene-divinylbenzene) Trimethylammonium
Particle Size:	10µm
Exchange Capacity:	0.19 ± 0.02meq/g
Mobile Phase Limits:	pH 2–12, 0–100% organic modifier

Anion/R Columns

Packing	i.d. x Length	PEEK Part No.	SS Part No.
Anion/R	4.1 x 100mm	—	269031
	4.1 x 250mm	—	269029
	4.6 x 100mm	269036	—
	4.6 x 150mm	269034	—

related product

Save money with IC All-Guard™ Kits on page 98. Kit includes one holder and three cartridges.

Guards for Anion/S and Anion/R Columns

Description	Qty.	Part No.
GA-1 Anion Cartridges	3	38108
All-Guard™ Cartridge Holder*	ea	80101

*Direct-Connect™ column coupler included.


Alltech® Cation Columns

Universal Cation and Universal Cation HR IC Columns

- Separate Groups I and II cations in one isocratic run
- Separate transition metals without post-column reaction
- New Universal Cation HR 3µm particle size for improved efficiency

Compatible with a variety of mobile phases including complexing acids (citric acid, tartaric acid, oxalic acid), mineral acids (nitric acid, hydrochloric acid, sulfuric acid) and non-complexing organic acids (methanesulfonic acid).

Alltech

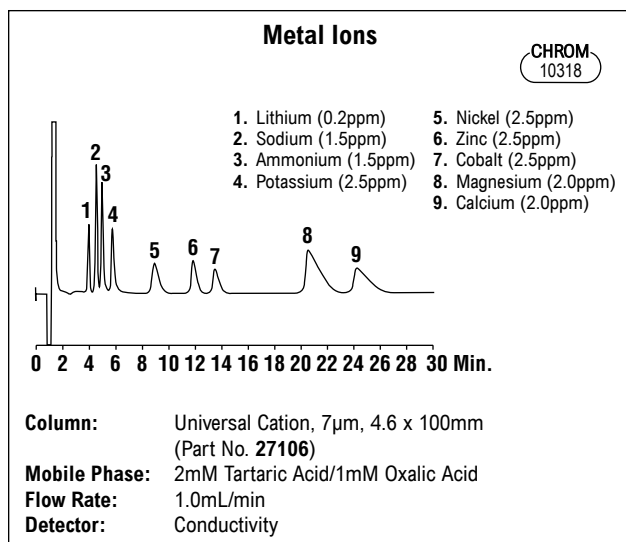

Universal Cation Specifications

Composition:	Silica coated with polybutadiene/maleic acid copolymer
Particle Size:	3µm and 7µm, spherical
Mobile Phase Limits:	pH 2–7, 0–100% organic modifier

Universal Cation Columns



Packing	i.d. x Length	PEEK Part No.	SS Part No.
7µm	4.6 x 100mm	27106	27100
3µm HR	4.6 x 100mm	—	23100
	7.0 x 53mm	—	23109

technical assistance

Contact Tech Support: Phone: 1.800.255.8324 (North America)
 Email: contact.alltech@grace.com
 Online: www.discoverysciences.com

more applications

To view our complete searchable chromatogram database visit www.discoverysciences.com/chromdb/

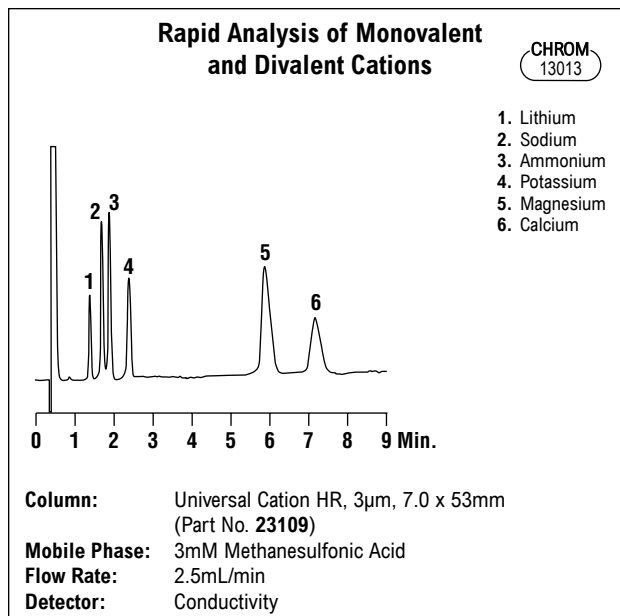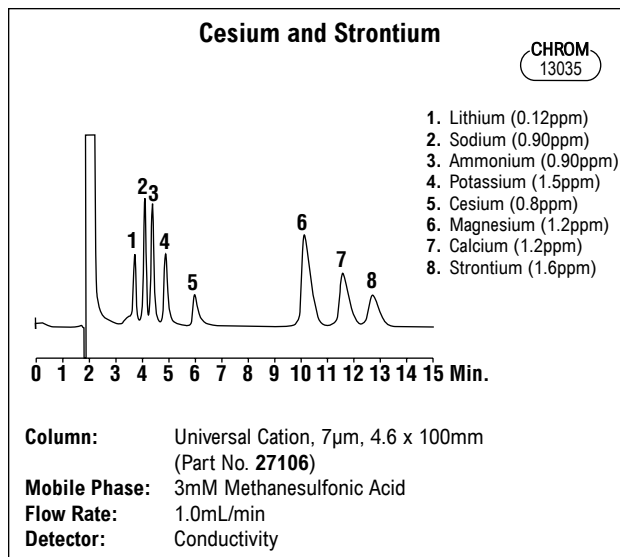

Guards for Universal Cation Columns

Description	Qty.	Part No.
GC-2 Universal Cation	3	271010
GC-3 Universal Cation HR	3	23110
All-Guard™ Cartridge Holder*	ea	80101

*Direct-Connect™ column coupler included.


related product

Save money with IC All-Guard™ Kits on page 98. Kit includes one holder and three cartridges.

Alltech® Anion Exclusion Columns

- Separate organic acids and weak acid anions
- Polymer-based for broad pH stability
- Economical choice for small inorganics

Anion Exclusion Columns separate organic acids and weakly ionized anions by an anion exclusion mechanism. Dilute mineral acids are the typical mobile phases for the separation. Acetonitrile may be used as an organic modifier to decrease the retention of hydrophobic compounds.

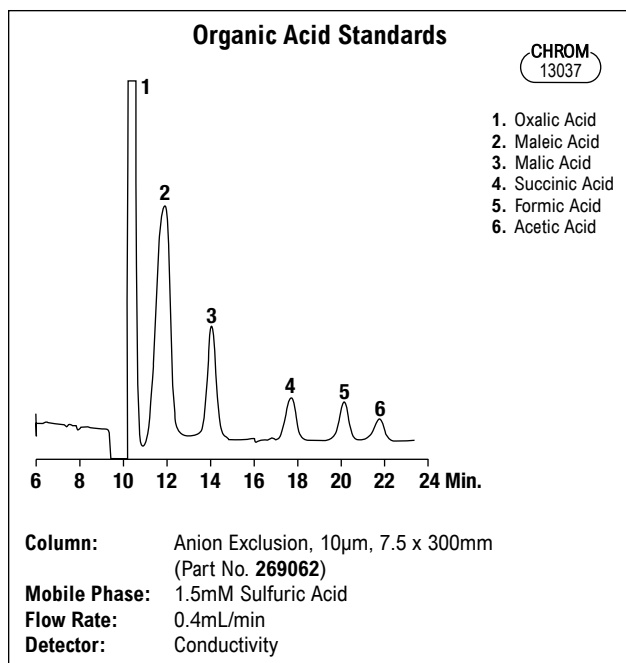

Anion Exclusion Specifications

Composition:	Highly sulfonated poly (styrene-divinylbenzene) cation exchanger
Particle Size:	10µm
Mobile Phase Limits:	<10% acetonitrile, no methanol <5% isopropyl or ethyl alcohols

Anion Exclusion Columns

Packing	i.d. x Length	SS Part No.
Anion Exclusion, 10µm	7.8 x 100mm	269068
	7.8 x 300mm	269006

Anion Exclusion All-Guard™ Cartridges

Description	Qty.	Part No.
Anion Exclusion GIE-1	3	38118
All-Guard™ Cartridge Holder*	ea	80101

*Direct-Connect™ column coupler included.

related product

Save money with IC All-Guard™ Kits on page 98. Kit includes one holder and three cartridges.


Alltech® Surfactant Columns

Surfactant/C IC Columns

- For non-suppressed analyses
- Silica-based for symmetrical peak shapes

Choose Surfactant C8 for short chain anionic surfactants such as alkyl sulfates and alkyl sulfonates.

Choose Surfactant C18 for long chain or aromatic surfactants such as xylene sulfonate and dodecyl benzene sulfonate.

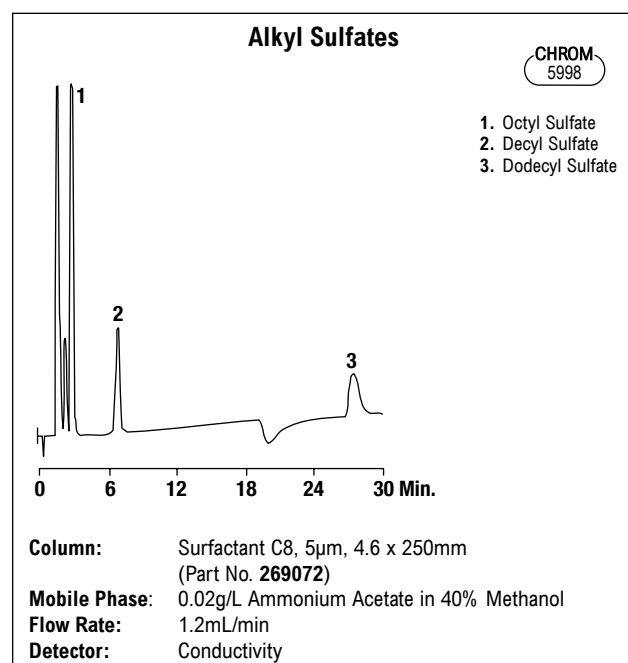

Surfactant/C Specifications

Composition:	Silica-based C8 or C18
Particle Size:	5 μ m
Mobile Phase Limits:	pH 2–7, 0–100% organic modifier

Surfactant/C Columns

Packing	i.d. x Length	SS Part No.
C8, 5 μ m	4.6 x 250mm	269092
C18, 5 μ m	4.6 x 250mm	269091

Guards for Surfactant/C Columns

Description	Qty.	Part No.
GRP-2 Revised Phase	3	38129
All-Guard™ Cartridge Holder*	ea	80101

*Direct-Connect™ column coupler included.