

Amides

Column: **MEGA-ACID (FFAP) - 0.53mm, 1.00µm, 10m**

Catalog Code: C-ACID-053-100-10

Conditions:

Injection: Split 250°C, 1.0µL.

Detector: FID 250°C.


Oven Temp.: 70°C, 10°C/min, 200°C.

Carrier Gas: Helium, 4 mL/min.

Peaks:

1. 3-picoline
2. N,N-dimethylformamide
3. N,N-dimethylacetamide
4. N-methylacetamide
5. Acetamide



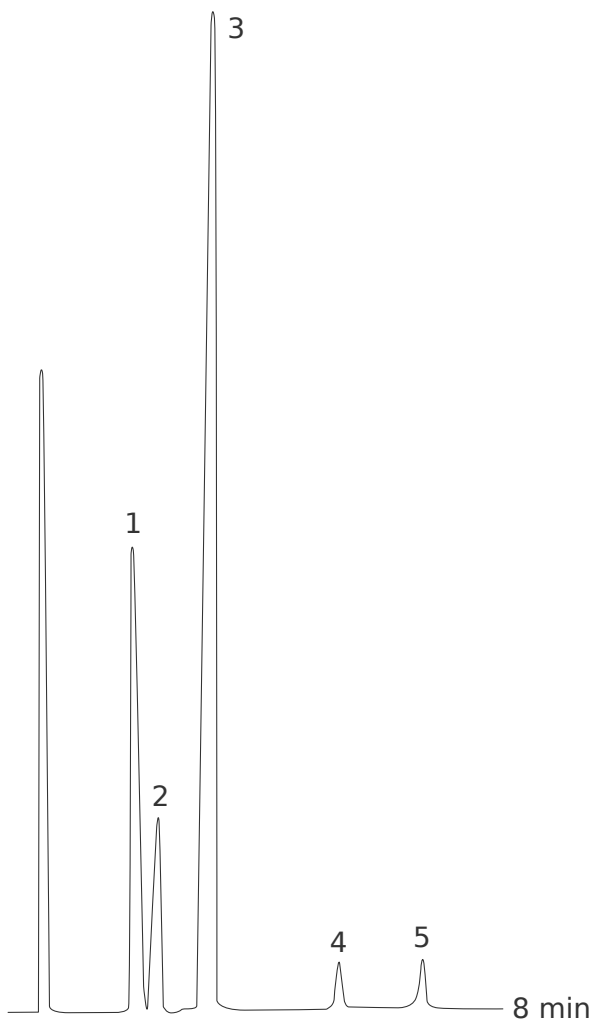

Amides

Column: **MEGA-ACID (FFAP) - 0.53mm, 1.00µm, 10m**

Catalog Code: C-ACID-053-100-10

Conditions:

Injection: Split 250°C, 1.0µL.

Detector: FID 250°C.

Oven Temp.: 70°C, 10°C/min, 200°C.

Carrier Gas: Helium, 4 mL/min.

Peaks:

1. 3-picoline
2. N,N-dimethylformamide
3. N,N-dimethylacetamide
4. N-methylacetamide
5. Acetamide

